


Optimal Exposure For Large Patients



- Designed for patients with a BMI 35 or above
- Unique elongated curved rack follows anatomic shape of the chest
- Extra deep blades open sternum with a gentle radius to minimize trauma

BMI 35 Large Adult Sternal Retractor with Interchangeable Atrial and Aortic Rakes

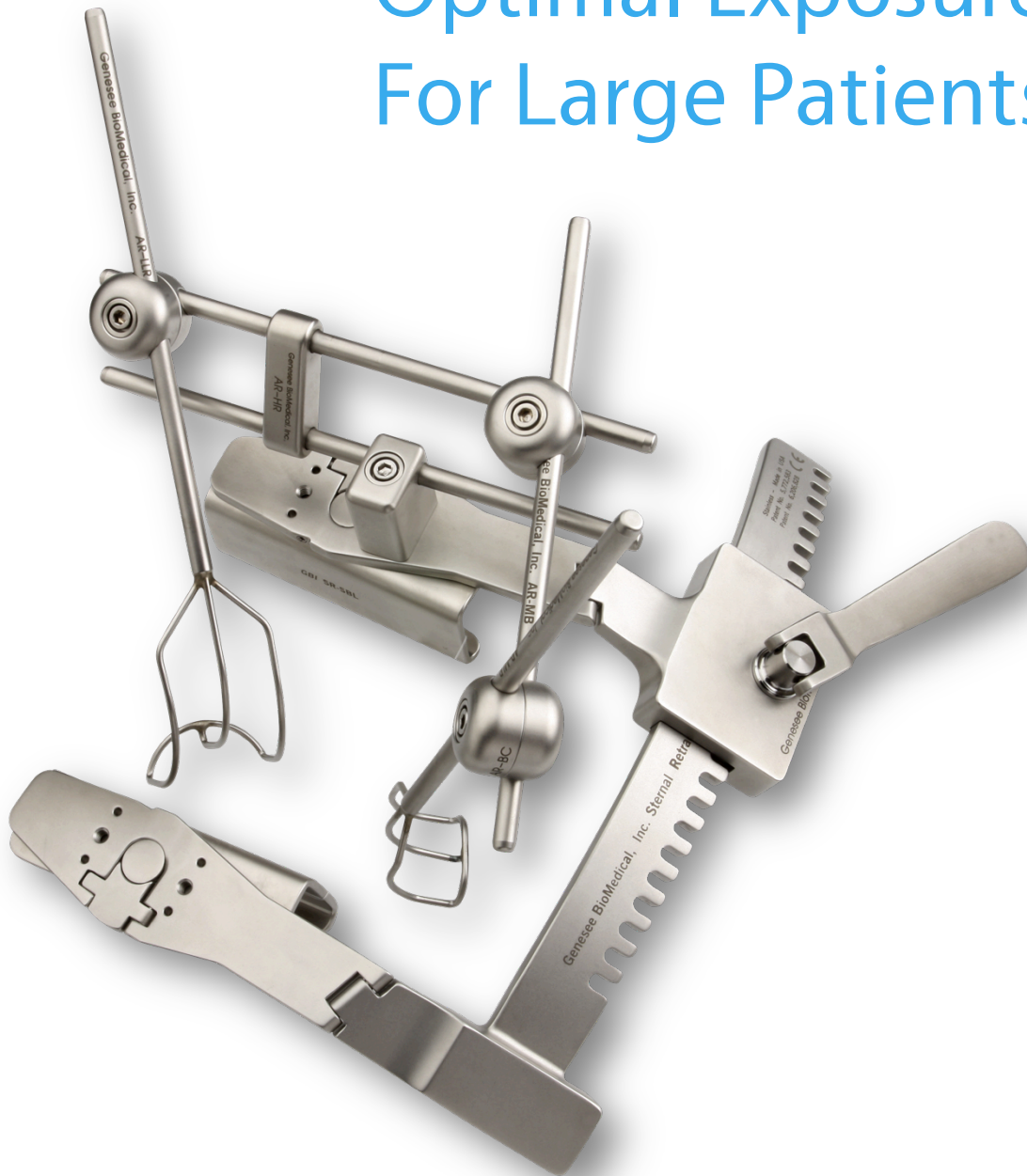

Unique Curved Rack Design

Innovative Design

- Unique rack design opens with a gently curved radius, optimizing exposure with minimal trauma to the chest and ribcage
- Anatomic curved shape distributes opening pressure equally on the sternum
- Pivoting blades equalize pressure points across the sternum
- Extra deep blades provide maximum exposure for patients with a BMI of 35 and above
- Suitable for adult patients above 90 kilos

Versatile System

- Unique ball clamps, H bar mounting rods and pillar clamps provide fast set up
- Interchangeable rakes provide excellent exposure and ease of use
- Changeable blades allow adjustment for individual patient

Ordering Information

For Customer Service; call toll free 800-786-4890 or 303-777-3000

Catalog Number

SRLP-AR

SRLP-1

ST-ST5

AAR-1

SOP

SR-DBL

SR-XDBL

Product Description

Sternal Latch Retractor with Atrial and Aortic Accessories

Large Patient Sternal Latch Retractor (Rack and Blades only) (Rack is 313 mm long, Arm is 191mm)

Standard Sterilization Tray

Atrial and Aortic Accessories

Suture Organizer (reusable)

Replacement Blades for SRLP

Adult Deep Latch Blade

(89mm wide by 38mm deep)

Adult Extra Deep Latch Blade

(89mm wide by 51 mm deep)



Genesee
BIOMEDICAL, INC.

Design beyond standard.

U.S. Patents #5,772,583 and #6,206,828.
Genesee BioMedical reserves the right to change or modify specifications without notice.
© 2015 Genesee BioMedical, Inc. All Rights Reserved. Printed in USA

www.geneseebiomedical.com